

Importing DICOM Data

Input Type: Database

The Database input option is for the import of DICOM image data into Analyze 14.0. The Database option enables the indexing of collections of DICOM images through a local database file. To create a database and for DICOM data import:

Open Input/Output.


Set Input Type to Database [1].

The first time the Database option is selected, a window will be displayed, asking:

‘There is no input database connection. What would you like to do?’

Select Create a new database. [2].


Note: If you have used previous versions of Analyze or AnalyzePro and have created DICOM databases you can use the ‘Browse filesystem for an input database’ to navigate to and select a .adb file. If this window is not returned or if this is not the first time you have opened the Database option, select Admin button [3] to open the Administer Databases window and create a new database.

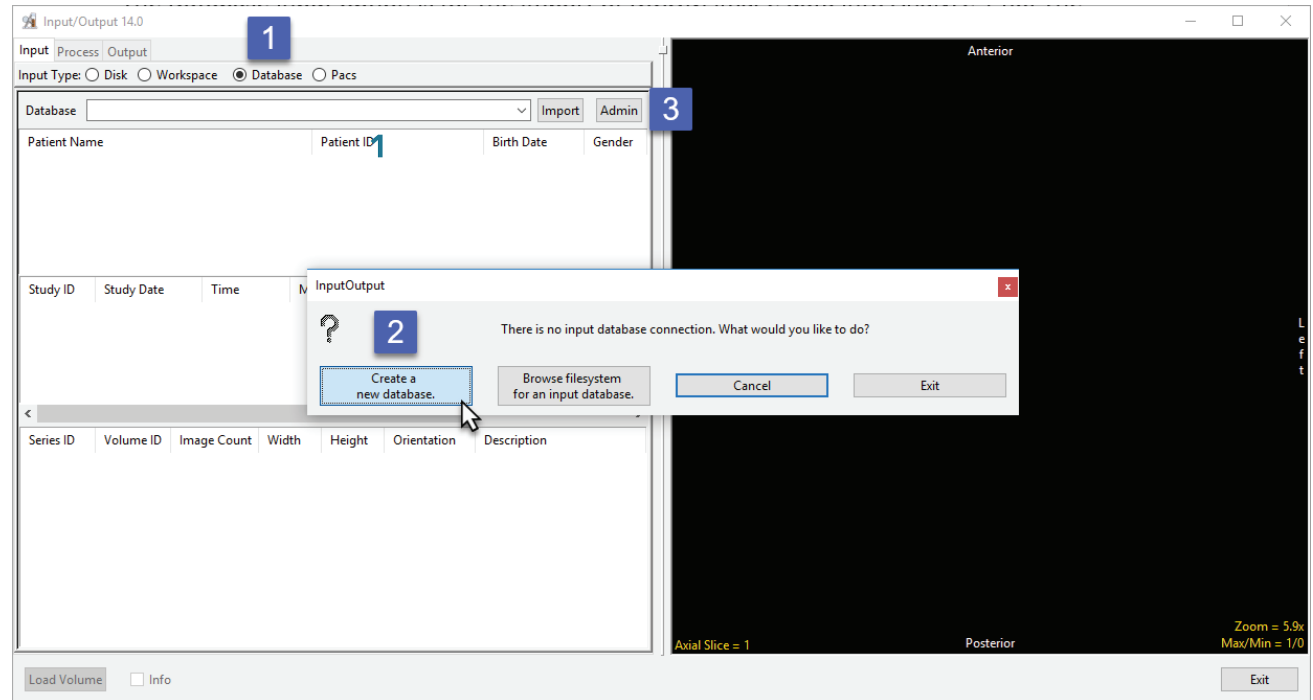

Importing DICOM Data (continued)

The Administer Database(s) window will open.

Name the database [3] and set the location where the database will be saved. [4].

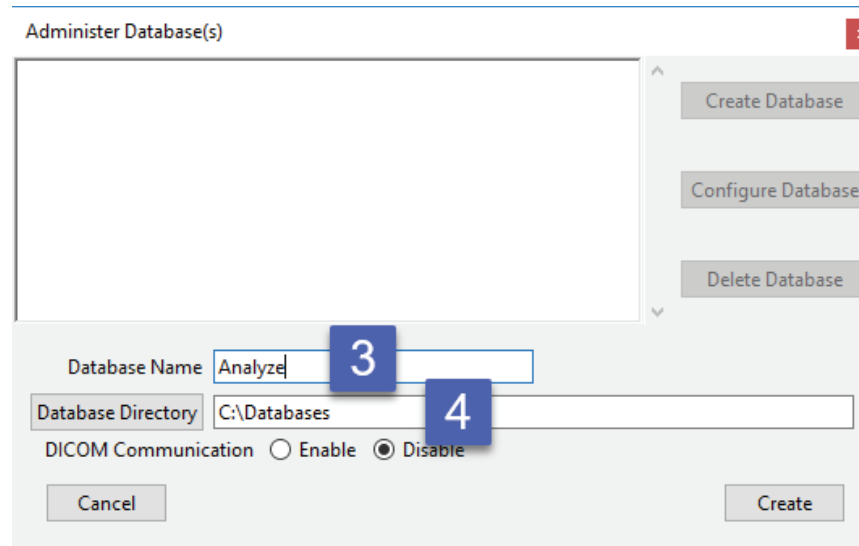

Importing DICOM Data (continued)

Enabling DICOM Communications (optional)

A DICOM receiver can be configured for the database. The receiver will run as a Windows service or Unix daemon and will actively listen on a specified port for DICOM image data to be pushed to it from PACS or any other DICOM node. To set up a DICOM receiver:

Set the DICOM Communication to Enable. [5].


Set the Port number [6] and Application Entity Title (AET) Name. [7]

Note that the AET Name should not go over 16 alphanumeric characters.

Make a note of the Port number and AET Name. Your PACS administrator will need this information along with the systems IP address to configure DICOM push from PACS.

Once configured on PACS you will be able to push DICOM image data from PACS to the Analyze 14.0 database. Note that if DICOM data is not being received after pushing from PACS you may need to open the port number for the Receiver on your system's firewall.

Click Create. [8].

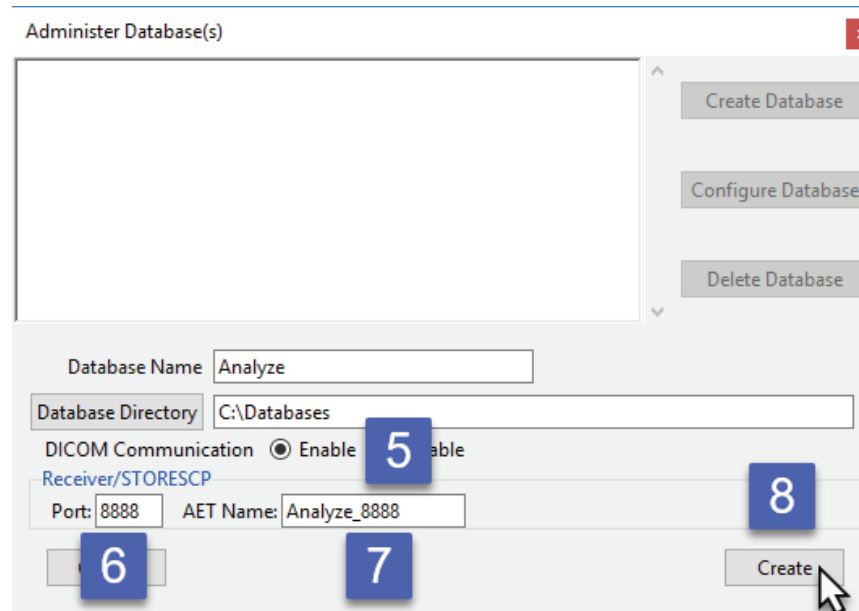


Importing DICOM Data (continued)

For Windows users an informational window will be returned informing your that a system prompt will follow to install services for the database. Click OK. [9]

In the Windows User Account Control prompt returned click Yes [10] to allow the program to install the receiver service.

Dismiss the service installed confirmation message by clicking OK [11].

Details about the receiver, as well as the associated database server, configured for the database are displayed at the bottom of the Administer Database(s) window [12]. Close the window.

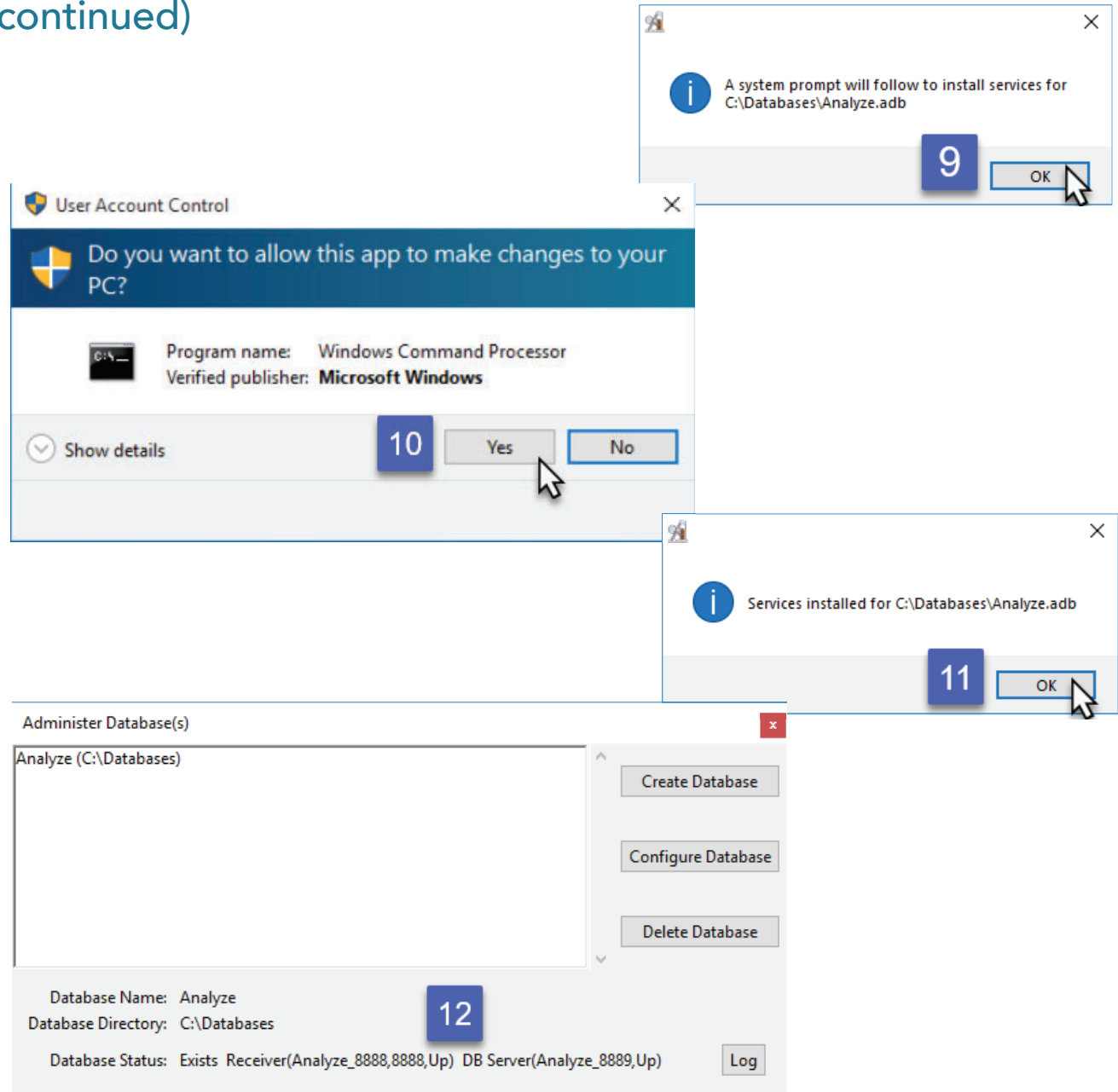

Importing DICOM Data (continued)




Importing Images

Once the Database is configured you can now start importing DICOM image data.

Click the Import button [13] or right-click in the window and select Import DICOM Images. [14]

Use the Directory Selector to navigate to the location of the DICOM images to import and click OK. [15]

A window will open stating Import all DICOM files found below <the specified folder>; click Yes. [16]

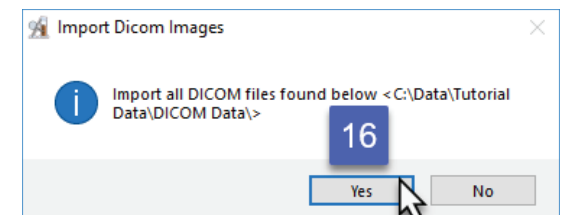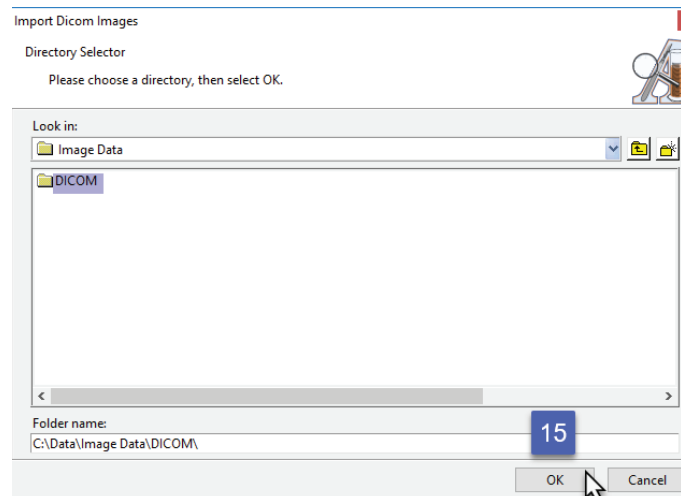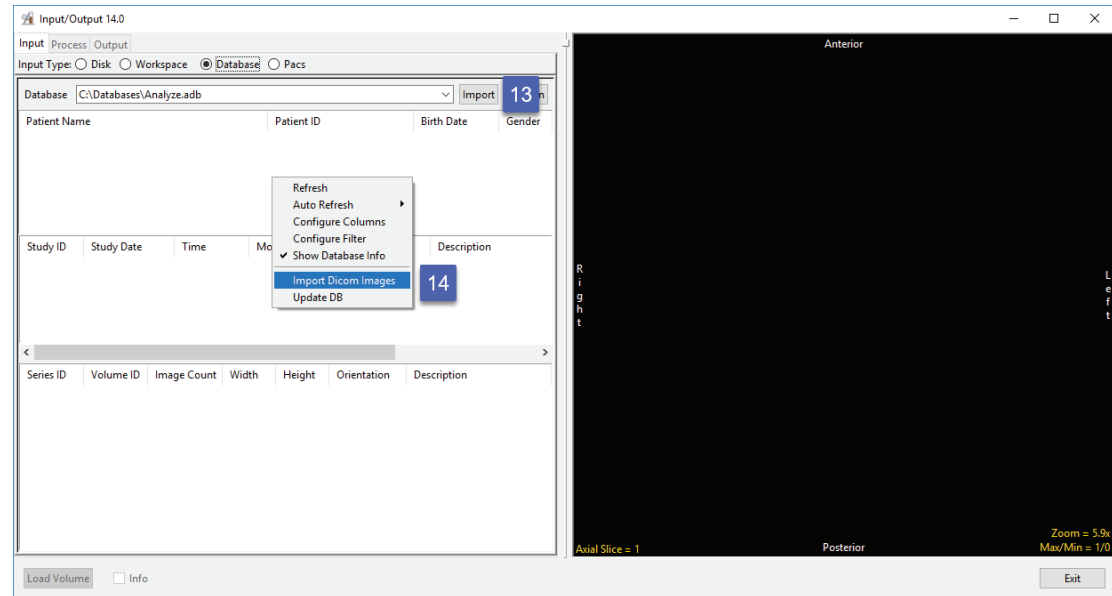

Importing DICOM Data (continued)

Loading DICOM Images

Once import is complete, the image data will be available in the database.

To load a data set into Analyze 14.0, select a data set and then click Load Volume. [17]

The screenshot shows the '11_Std11_Ser200_V1 - Input/Output' window. The 'Input' tab is active, and 'Database' is selected as the input type. The database path is 'C:/Databases/Analyze.adb'. Below this is a table of patient data:

Patient Name	Patient ID	Birth Date	Gender
MR SIGNA CONTOUR 0.5T	11	NA	F
MR SIGNA LX 1.0T	337	NA	F
MR SIGNA LX 1.0T	930	NA	F
MR SIGNA LX 1.5T	28	NA	O
MR SIGNA PROFILE 0.2T	93	NA	M
MR dataset Brain	reg2	Jan 01, 1911	M
ULTRASOUND GE	000-00-0000	NA	F

Below the patient table is a table of study data:

Study ID	Study Date	Time	Modality	Referring Phys	Description
11	Apr 02, 1997	13:57:45	MR	MELOFF	e+1 MR BRAIN

At the bottom of the window, there is a table of series data:

Series ID	Volume ID	Image Count	Width	Height	Orientation	Description
200	1	18	256	256	Axial	C5,1,4,5,0,7,1,1,0,85,108
200	2	18	256	256	Axial	C5,1,4,5,0,7,1,1,0,85,108

To the right of the tables is a large window displaying an axial MRI brain scan. The scan is labeled 'Anterior' at the top and 'Posterior' at the bottom. The left side is labeled 'Right' and the right side is labeled 'Left'. The scan shows a cross-section of the brain with a zoom level of 2.3x and an axial slice of 9. The maximum and minimum values are 993/0. At the bottom left of the scan window, there is a 'Load Volume' button with a mouse cursor over it, and an 'Info' checkbox. At the bottom right of the scan window, there is an 'Exit' button. A blue box with the number '17' is overlaid on the 'Load Volume' button.