


Storing Data in a Database

- Select the Database option [1].
- If required, select a database from the drop-down [2].
- Click Store Volume [3].

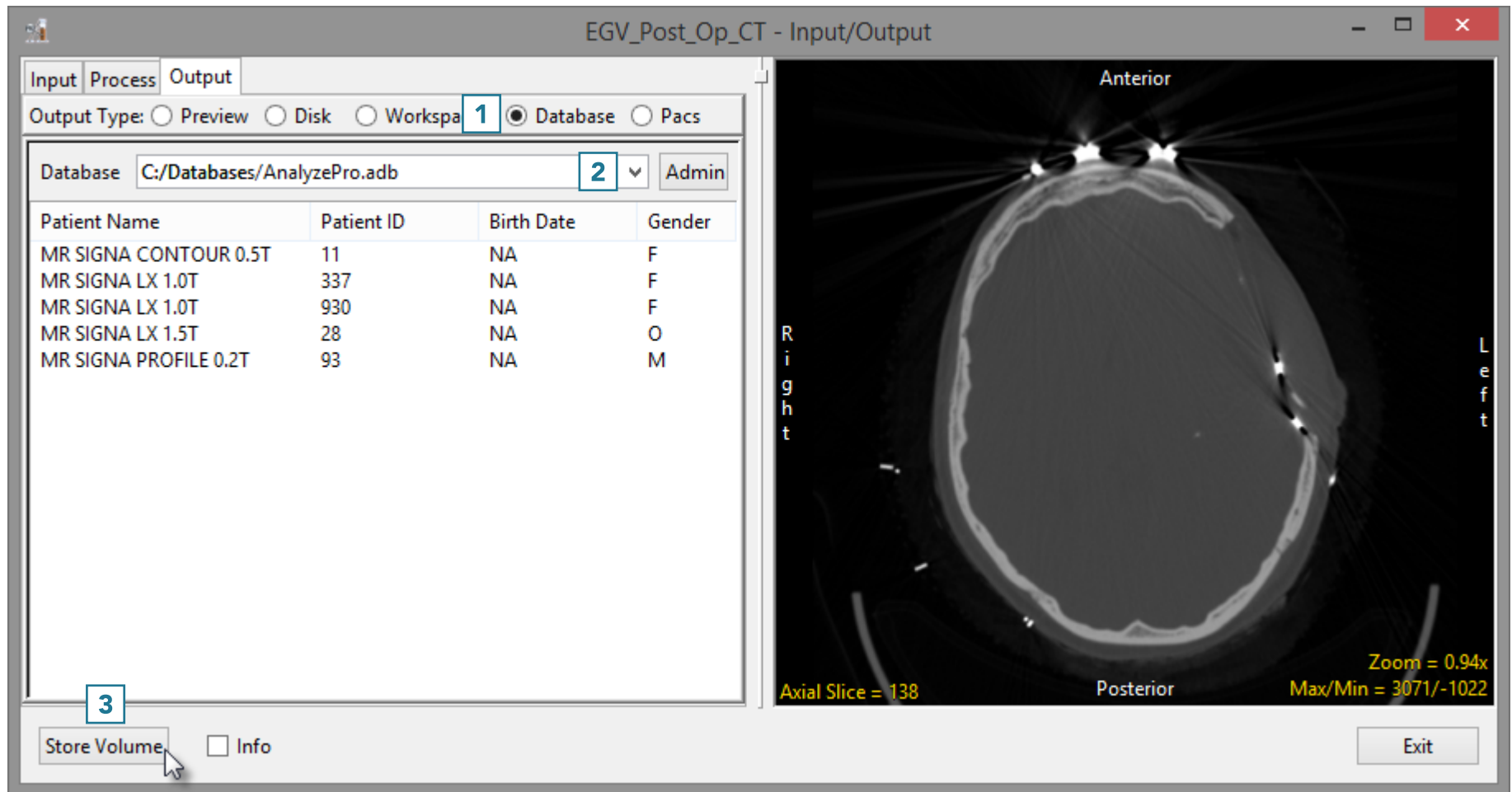

Storing Data in a Database (continued)

- The Edit/Add DICOM tags window will open. Use this tool to add or edit any DICOM tag information.
- Click Apply Tags and Continue [4].
- The data will be exported to and merged into the database.

Edit/Add Dicom tags

Edit DICOM tags

Patient

ID (0010,0020) 500

Name (0010,0010) Joe Bloggs

DOB (0010,0030)

Gender (0010,0040) M

Study

ID (0020,0010)

Description (0008,1030)

Date (0008,0020)

Accession (0008,0050)

Modality (0008,0060)

Portion (0008,0015)

(0020,000d) Make New Study Instance UID

Series

ID (0020,0011)

Orientation Axial Coronal Sagittal

Description (0008,103e)

Procname (0008,1030)

(0020,000e) Make New Series Instance UID

Other Tags

(gggg,eeee) |

4

Apply Tags and Continue Clear Add Tag Cancel